


There are two ways to fix the matter which regarding to the new format of EDF file


1. Run the tool EdfHeaderConverter.exe and select the files which need to convert. This way will be run for eachtime when you want to run a new file
 2. Copy getfiletype.m and sopen.m files to 'plugins\Biosig3.1.0\biosig\t200_FileAccess'. With this way you can run edf file with old or new Emotiv EDF header format.
- At first, the customer need to open Matlab and run eeglab.

In EEGLab, click File > Import data > Using EEGLab functions and plugin > From EDF/EDF+/GDF files...

Then will show a message to install Biosig. Please click Yes to install

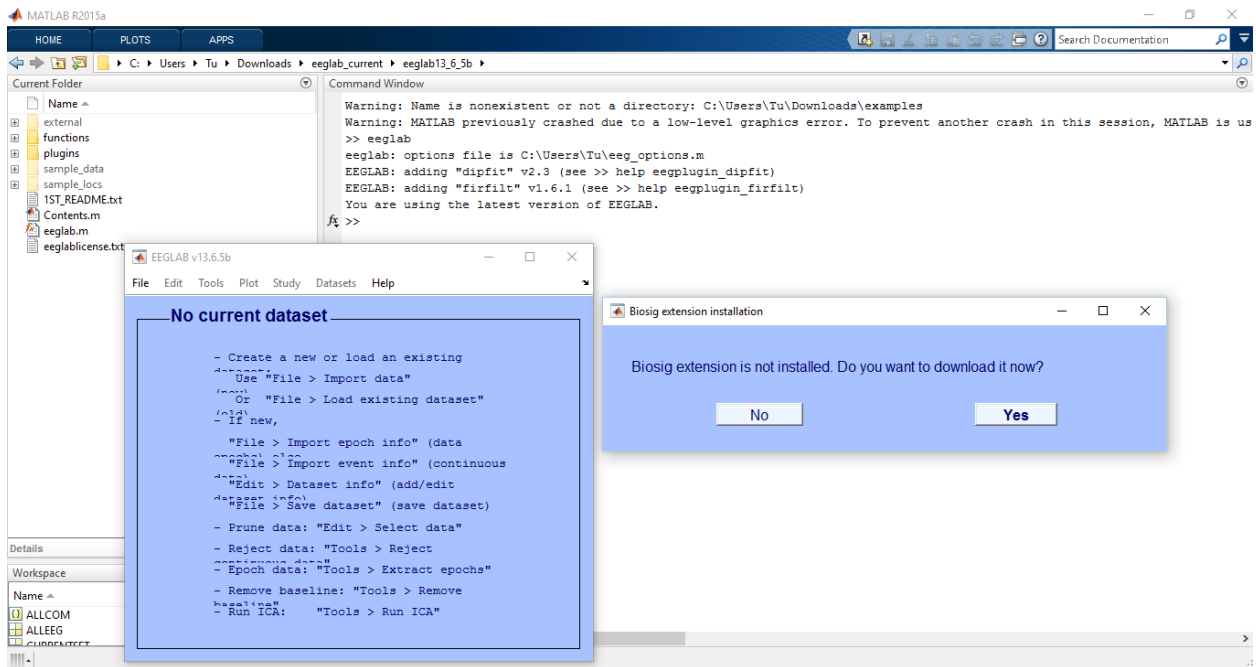

Click Yes if there is a pop up warning appear

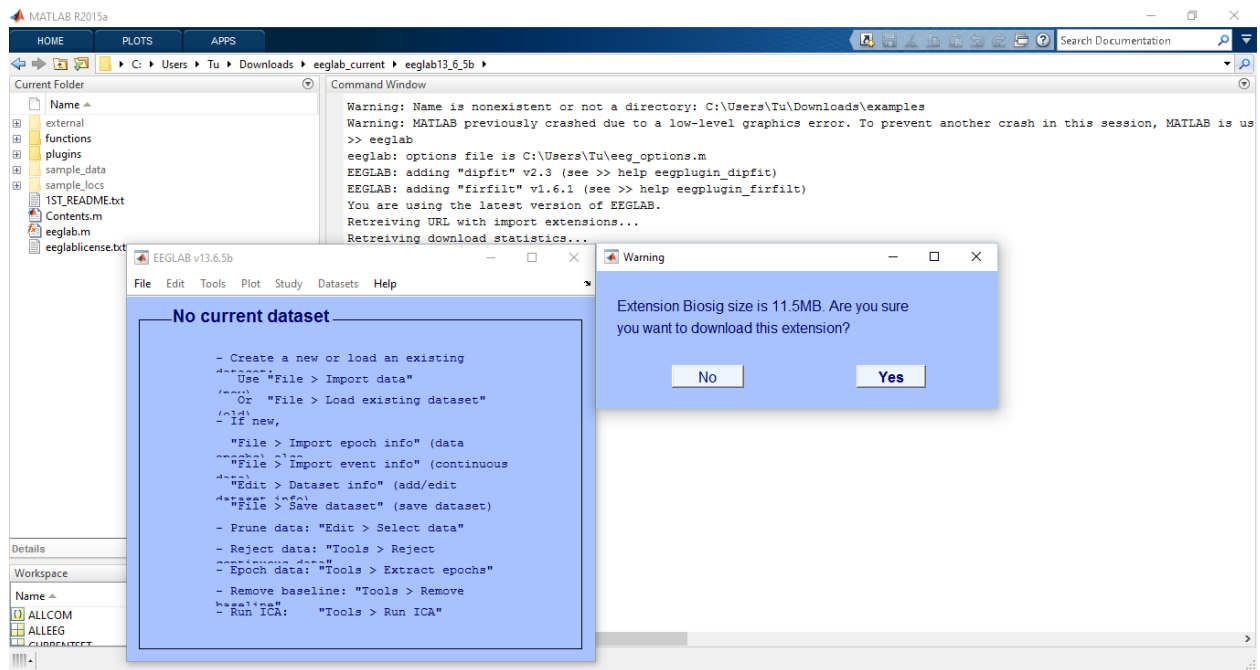

- Copy getfiletype.m and sopen.m files to 'plugins\Biosig3.1.0\biosig\t200_FileAccess'
- Then you can select file to add for eeglab